



**EKOTEKS LABORATUVAR ve GÖZETİM
HİZMETLERİ A.Ş.**
Esenyurt Firuzköy Bulvarı No:29 34325 Avcılar
İstanbul/ TÜRKİYE



TEST REPORT
DENEY RAPORU


AB-0583-T
20024672
07-20

Customer name: PLATİN MODA TASARIM TEKS. SAN. VE DIŞ TİC. LTD. ŞTİ.
Address: Çobançeşme Mah. Selvi Sokak No:3/K:3 BAHÇELİEVLER/ İSTANBUL
Buyer name: PLATİN MODA TASARIM TEKS. SAN. VE DIŞ TİC. LTD. ŞTİ.
Contact Person: ABDULLAH YİĞİT
Order No: -
Article No: -
Name and identity of test item: Blue surgical un-woven gown . (Claimed to be; DISPOSABLE GOWN NON STERILE-LEVEL 3 ONE SIZE FIT MOST 100% POLYPROPYLENE)
The date of receipt of test item: 17.07.2020
Re-submitted/re-confirmation date: -
Date of test: 17.07.2020-17.07.2020
Remarks: -
Sampling: The results given in this report belong to the received sample by vendor.
End-Use: -
Care Label: Not specified.
Number of pages of the report: 3

The Turkish Accreditation Agency (TURKAK) is signatory to the multilateral agreements of the European co-operation for the Accreditation (EA) and of the International Laboratory Accreditation (ILAC) for the Mutual recognition of test reports.

EKOTEKS LABORATUVAR ve GÖZETİM HİZMETLERİ A.Ş. accredited by TÜRKAK under registration number [AB-0583-T] for ISO 17025:2017 as test laboratory.

The test and/or measurement results, the uncertainties (if applicable) with confidence probability and test methods are given on the following pages which are part of this report.

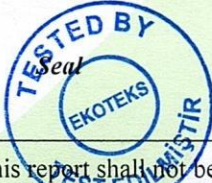

Date
18.07.2020

Customer Representative
Özlem ULUS

Head of Testing Laboratory
Sevim A. RAZAK
18.07.2020

This report shall not be reproduced other than in full except with the permission of the laboratory.
Testing reports without signature and seal are not valid.

**EKOTEKS LABORATUVAR ve GÖZETİM
HİZMETLERİ A.Ş.**

AB-0583-T

20024672

07-20

REQUIRED TESTS	RESULT	COMMENTS
PHYSICAL PROPERTIES		
Water Permeability	P	
Water Resistance (Impact Penetration)	P	
P: Pass F: Fail R: Refer to retailer technologist. Evaluation was performed according to the AAMIPB70 guidelines.		

REMARK: Original samples are kept for 3 months and all technical records are kept for 5 years unless otherwise specified. If requested, measurement uncertainty will be reported. But unless otherwise specified, measurement uncertainty is not considered while stating compliance with specification or limit values. The reported uncertainty is based on a standard uncertainty multiplied by a coverage factor $k=2$, providing a level of confidence of approximately 95 %. Tests marked (*) in this report are not included in the accreditation schedule.



Gen.f136-2/03

This report shall not be reproduced other than in full except with the permission of the laboratory.
Testing reports without signature and seal are not valid.

**EKOTEKS LABORATUVAR ve GÖZETİM
HİZMETLERİ A.Ş.**

AB-0583-T

20024672

07-20

TEST REPORTS

WATER PERMEABILITY; AATCC 127:2017-OPTION 2

Hydrostatic Head Tester, Textest marka Fx 3000 model
Temperature of water 20°C. Pressure increase ratio 10 mbar/dk.
Performed in the conditioned room (21±2°C-65%±5)

	<u>RESULT</u>	<u>REQUIREMENT</u>
Sample 1	213,1cmSS	
Sample 2	197,8 cmSS	
Sample 3	198,9 cmSS	
Average	203,2 cmSS	≥ 50 cmSS (Level 3)

WATER RESISTANCE-IMPACT PENETRATION TEST; AATCC 42:2017(*)

Test samples were conditioned for at least 4 hours prior to testing (65 ± 5% RH, 21 ± 2 ° C).
Water temperature : 27±1°C
Amount of sprayed water: 500 ml

	<u>RESULT</u>	<u>REQUIREMENT</u>
Water penetration (g)	0.0 g	≤1,0 g (Level 3)